


PIGNAT

# Continuous distillation

Référence : DVI/100

## Study of the continuous distillation

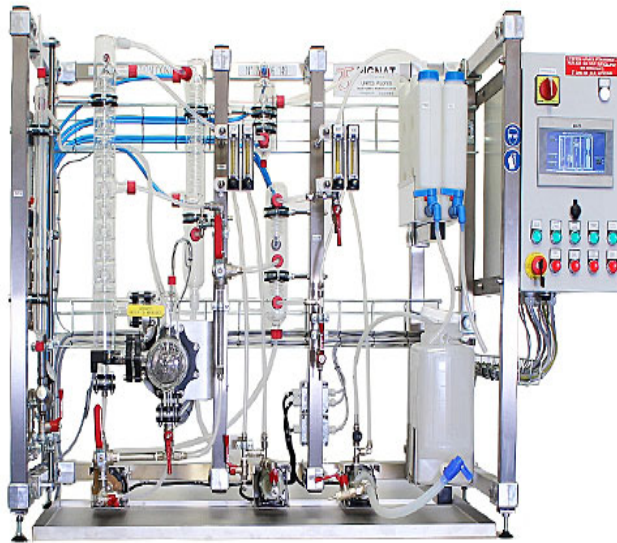

Dim: 165 x 70 x 115 cm – 150 kg

Thermosiphon boiler

Bubble trays column

External reflux

### GENERAL SPECIFICATIONS

- Feeding pump, with by- pass valve.
- Distillate pump.
- Column bottom pump.
- Glass preheating, two coils.
- Glass boiler, heater.
- Low level detector.
- Glass column bubble cap plates.
- Glass condenser, glass coil.
- Cooler, in glass, glass coil.
- Distillate receiver.
- Bottoms receiver.

### Instrumentation

- Temperature probes Pt100Ω.
- Differential pressure.
- Flowmeters.
- Level sensor with vibrating blades.



**REPRESENTED BY:**  
Allegheny Educational Systems, Inc.  
320 East 3rd Avenue  
Tarentum, PA 15084  
Phone: 800-232-7600  
www.alleghenyedusys.com

# DESCRIPTION

---

## Distillation

Distillation is a process including heat and material transfers in order

to separate one or several products of relative different volatility.

The separation is carried out by a chaining of vaporizations and condensations.

The theoretical base is that a vapor produced from a mixture is in equilibrium with the mixture liquid phase and has generally a higher composition in the most volatile compound.

The separation is carried out in a column (packed or plates column) designed to get the best contact and material exchanges between the climbing vapor and the falling liquid into the column.

This separations can be done in continuous or discontinuous (batch) process.

## Study of column hydrodynamics

Determination of overflow conditions

## Study of main parameters

Mass Balance

Thermal Balance

Yield and hold-up

*A automated version is available,  
contact-us for more information*



Unités livrés avec un manuel pédagogique  
et dossier technique



Possibilité d'adapter les unités à vos  
besoins



Mise en service, formation des  
formateurs